




Magnetic Pulsar^{MP6}

Troubleshooting



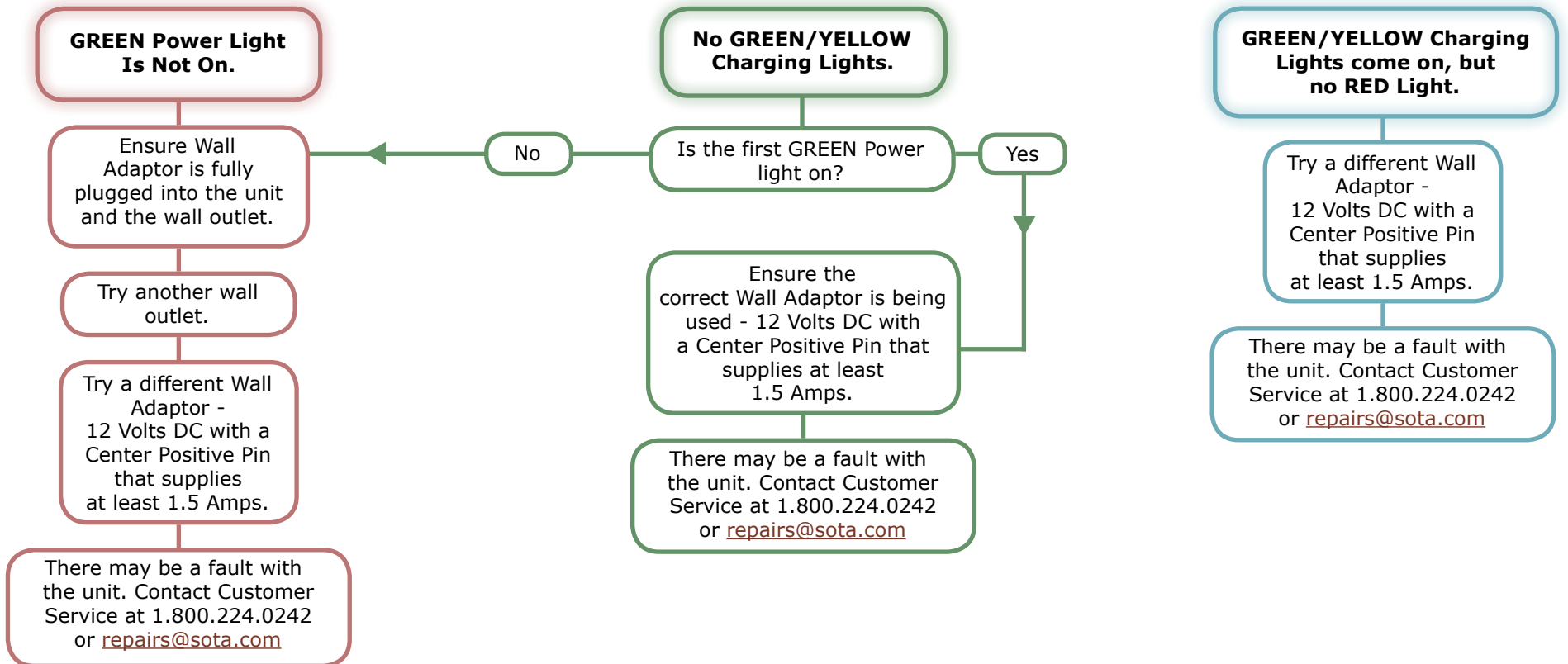

Magnetic Pulsar^{MP6}

Troubleshooting

Pulse Takes Too Long To Build and Release

Ensure the correct Wall Adaptor is being used - 12 Volts DC with a Center Positive Pin that supplies at least 1.5 Amps.

Try a different Wall Adaptor - 12 Volts DC with a Center Positive Pin that supplies at least 1.5 Amps.

There may be a fault with the unit. Contact Customer Service at 1.800.224.0242 or repairs@sota.com

Hand Paddle heats up after one or two cycles.

This is normal. The Hand Paddle contains a coil through which the magnetic field is released causing the coil to generate heat.

If the heat from the paddle is uncomfortable we suggest:

Put a cloth around the paddle - a sock works well.

Turn the unit off and allow the paddle to cool before continuing.

Unit is beeping and the RED and GREEN/YELLOW Lights are flashing rapidly.

There is a fault with the unit. Contact Customer Service at 1.800.224.0242 or repairs@sota.com